

Kestrel 1000

Digital Scientific Frame Transfer EMCCD

128 x 128 • 24 μ m x 24 μ m pixels • Cooled to -20°C • 1800fps •



Key Features and Benefits

Ultra low noise readout with THE FASTEST speeds

- **Up to 1800 fps with ROI, 1000 fps in full resolution**

The fastest EMCCD on the market

- **128 x 128 Back-thinned EMCCD sensor**

Enables optimum image resolution in low light imaging applications

- **16 bit Camera Link output**

Realtime imaging for low latency photon to digital image

- **Up to 95% QE from back-illuminated sensor**

Optimum Photon collection

- **Strong VIS and NIR response and ultrawide bandwidth**

From 350nm through to 1100nm

Resolution	128 x 128
Pixel Size	24μm x 24μm
Readout Noise	<1e
Frame Rate	1,000fps
Camera Link	16bit

Specification for Kestrel 1000

Sensor Type	1/8" Back Thinned Frame Transfer EMCCD
Active Pixel	128 x 128
Pixel Size	24µm x 24µm
Active Area	3.1mm x 3.1mm
Full Well Capacity	100ke-
Non-linearity	<1%
Readout Noise (RMS)	EM Gain ON: <1e- EM Gain OFF: 290e- Typical
Full Resolution Frame Rate	1,000fps
Dark Current (e/p/s)	<1 @ -20°C
Digital Output Format	16 bit Camera Link (base configuration)
Peak Quantum Efficiency	95%
Spectral Response	350 - 1100nm
TE Cooling	-20°C With Fan / liquid cooled @ 25°C Ambient
Binning	1x1 up to 2x2
Lens Mount	C-Mount
Synchronisation	Trigger IN and OUT - TTL compatible
Power Supply	12V DC ±10%
Total Power Consumption	<80W
Operating Case Temperature ¹	-20°C to +55°C
Storage Temperature	-30°C to +60°C
Dimensions (L*W*H) ²	121mm x 140mm x 113mm
Weight (no lens)	< 1.5kg

Raptor Photonics Limited reserves the right to change this document at any time without notice and disclaims liability for editorial, pictorial or typographical errors.

Ordering Information

Camera

Kestrel 1000 EMCCD digital camera KE60V1000-BV-CL

Power Supply Unit RPL-HR4-K

Optional Accessories

Mini PC with XCAP STD and frame grabber RPL-PC-mf2280

Thunderbolt frame grabber RPL-mf2280

EPIX® EB1 base CL card RPL-EPIX-EB1

EPIX® XCAP STD software RPL-XCAP-STD

MDR-SDR Camera Link Cable, 2m³ RPL-MCL-CBL-2M

Thermoelectric Water Chiller Unit RPL-CHILLER

Chiller Tubing RPL-WTUBE-NINOX

Optical Lenses⁴ RPL-xx-xxxx

Note 1: Extended operating temperature range on request.

Note 2: Dimensions include all connector parts on camera interface


Note 3: Longer CL cable available up to 25M

Note 4: Please consult us to check our range of lenses

Demo is available on request.
Pricing AOR subject to volumes.

Detailed technical drawings
can be downloaded at
www.raptorphotonics.com

Quantum Efficiency

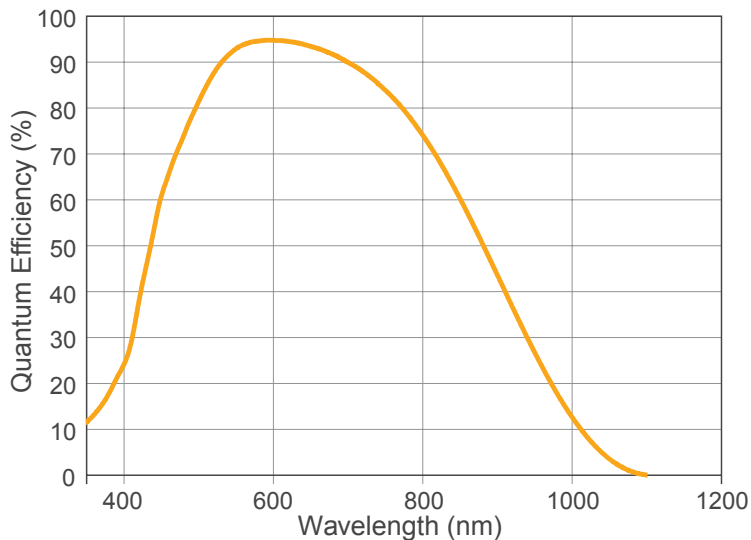

*Data supplied by sensor manufacturer

Applications

- Adaptive Optics and Astronomy
- Calcium signaling
- High resolution fluorescence imaging
- Hyperspectral imaging
- X-ray & High energy
- Particle Image Velocimetry
- High speed object tracking

Document #: USKE60V1000-BV-CL 0322